

Lab Equipment





Heathrow Scientific Equipment Family



MagFuge®

Did you know?

Heathrow Scientific has strict standards that include 3rd party plant reviews and 100% mutilevel product inspections. These actions result in a 99.7% reliability on our equipment and 99.99% reliability on our laboratory supplies.

Safety certified and approved Meets international standards.



High-speed Centrifuge is easy to use and designed to deliver maximum speed (12,500rpm/9,783x g) and flexibility with rapid acceleration and fast controlled breaking

Ideal for most protocols requiring fast spins including PCR, microfilter cell separation and HPLC. Versatility to handle 0.2, 0.5, 1.5, 2.0 and 5.0 mL tube sizes

Intuitive One-Dial Programming

- Set speed and time with ease:
- Adjustable speed from 500 to 12,500 rpm
- Speed setting for RPM or RCF
- Adjustable time from 30 seconds to 60 minutes

Safe to Use

- Safety features help ensure it will not operate when the lid is open, a rotor is jammed, or the unit is out of balance

Item No.	Color	L x W x H in	L x W x H cm
120581	Grey/Purple 🔳 🔳	- 9.2 x 7.3 x 4.7	23.5 x 18.5 x 12
120582	Grey/Blue	9.2 X 7.3 X 4.7	23.3 × 10.3 × 12
Replaceme	nt parts		
120622	MagFuge [®] Rotor, 12 x 1.5/2.0 mL tube		
120623	MagFuge® Rotor, 6 x 5.0 mL tube		
120624	MagFuge® Stirrer Rotor		
120682	MagFuge® Silicone Mat		
120527	Power Adapter		
120015	Tube Adapters, 0.2 mL, 0.5 mL, Pack of 12 (6 ea. 2 Packs required)		
HS100502	Universal Power Cord Set		

Patent Information

Community Design No. 003515337-0001-0002 Euro. Pat. App. 3 246 088 US Design Patent No. D814,045 US Patent No. 10471439



High-performance Magnetic Stirrer is designed using high quality rare earth magnets creating a robust magnetic coupling that significantly reduce the chance of spin-outs

Perfect for work with sensitive samples and protocols up to 3L requiring aggressive mixing for full dilution

Easy to Use

- Set stir mode, speed and time with easy one-dial programming:
- Select Mode: clockwise, c. clockwise and oscillate
- Adjustable speed from 50 to 2,500 rpm
- Continuous and adjustable time settings from 30 seconds to 60 minutes

Durable Accessories Included

 Re-usable, easy to clean silicone mats included to keep vessel stable and prevent spills and two NdFeB stir bars (40 x 12 mm) with chemical resistant PTFE coating included

Package Includes		
Centrifuge/Stirrer	Stirrer Rotor	12 Tube Adapters 0.2 mL & 0.5 mL
0		· · · · · · · · · · · · · · · · · · ·
5.0 mL Tube Rotor	2 Stir Bars, 35 m	m x 12 mm, cylindrical shape
1.5/2.0 mL Tube Ro	tor 2 Silicone Mats f	or Stirrer Lid
1 Low Voltage, Doub	ole Insulated Power Adapt	er (100-240 VAC, 50/60 Hz) with 4
Interchangeable Plu		
Specifications: Cer	ntrifuge	
Speed Range	500 – 999 rpm (increme 5,000 – 12,500 rpm (inc	nts of 1); 1,000 – 5,000 rpm (increments of 5); rements of 25)
RCF	16 – 9,783 x g (max RCF	calculated)
Run Time Ranges 30 seconds - 60 minut		s (5 second increments)
Creations: Ma	anatia Ctinnan	
Specifications: Ma	gneuc Surrer	
Speed Range	50 – 999 rpm (increments of 1); 1,000 – 2,500 rpm (increments of 5)	
Run Time Range	Stirrer Mode	30 seconds to 60 minutes or continuous (5 second increments)
	Oscillating Stirrer Mode	1 second – 60 minutes or continuous
	Oscillating Period	30 seconds – 15 minutes
Capacity	3 L plus, Diameter of ves	ssel not to exceed 130 mm

Quick-spin Button

Gusto® High-Speed Mini Centrifuge

Maximum speed with minimal space

This compact centrifuge is ideal for protocols calling for fast spins with a higher RCF requirement

Made of durable ABS with UV-Resistant coating

User Friendly

- Intuitive Push button programming and bright LCD screen make it easy to use

Quiet Operation

- Brush-less motor enables quiet operation and low vibration

Safety Features

- Lid safety feature prevents opening while rotor is turning and clear polycarbonate allows user to see rotor movement for added visual safety
- Includes automatic shut-off feature for sample imbalance and tilt

ltem No.	Color	L x W x H in	L x W x H cm
HS10050	Blue 🗖	8 x 6.7 x 4.5	20.3 x 17.1 x 11.4
Replaceme	nt parts		
HS100502	Universal Power Cord Set		
HS100503	Gusto® Rotor 12 x 1.5/2.0 mL Tube (with Cover and Knob)		
HS100504	Gusto® Rotor 4 x 8 Strip 0.2 mL Tube (with Knob)		
120015	Tube Adapters, 0.2, 0.5 mL, Pack of 12 (6 ea. 2 Packs required)		
HS100506	Gusto® Rotor Cover (for 12 x 1.5/2.0 mL Tube Rotor Only)		
120527	Power Adapter		

for momentary spin requirements,

Adjustable time and speed

Year Warranty

press and hold to operate, release to end spin cycle Back-lit Display showing time, speed and operational information on the same screen



Specifications	
Speed Range	500–12,500 rpm
Run Time Ranges	1 to 30 minutes
RCF	15 - 9,783 x g (max RCF calculated)



Package Includes
Gusto® High-Speed Mini Centrifuge
12-Place 1.5/2.0 mL Tube Rotor
32-Place 0.2 mL Tube Rotor (singles, or 4 x 8 Strip)
12 Tube Adapters 0.2 mL
12 Tube Adapters 0.5 mL
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs



Sprout[®] Plus Mini Centrifuge

Mini but Mighty - ideal for small runs

Easy to Use

- Tool-free SnapSpin rotors and quick-spin options make the centrifuges fast and easy to use

Ergonomic Design

- The lid was functionally designed to fit in the _ palm of the hand for one handed lock and unlock allowing for quick tube extraction
- Large lid size accommodates longer and skirted tubes and tall tube lids
- The housing is made from ABS plastic and lid is polycarbonate plastic enabling the unit to be easily cleaned

Safety Measures

Sprout[®] Plus Item No.

Replacement parts HS100011

120610

120611

HS100012

120148

120491

120015

Fitted with a braking function, the unit quickly stops spinning when powered off or when lid is opened _

A weighted, rubberized base eliminates potential movement and dampens vibration and noise







Package Includes

Specifications		
Speed Range	Approx. 6,000 rpm (fully loaded)	
RCF	2,000 x g	

Sprout® Plus Mini Centrifuge
6-Place 1.5/2.0 mL Tube Rotor
16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip)
6 Tube Adapters 0.2 mL
6 Tube Adapters 0.5 mL
4 Lid O-Rings
Storage Case for Rotors and Adapters
Tube Rack
1 Low Voltage, Double Insulated Power Adapter (100- 240 VAC, 50/60 Hz) with 4 Interchangeable Plugs

			10-Flace 0.2 IIIL Tube holdi (Sillyles, u	
Color	LxWxH		6 Tube Adapters 0.2 mL	
	in	cm	6 Tube Adapters 0.5 mL	
Green	——6 x 5 x 4.3	15.3 x 12.8 x 10.9	4 Lid O-Rings	
Blue	0 7 0 7 1.0	13.3 X 12.0 X 10.3		
ts			Storage Case for Rotors and Adapters	
Sprout [®] Plus Rotor SnapSpin 6 x 1	.5/2.0 mL Tube		Tube Rack	
Sprout® Plus Rotor SnapSpin 2 x 8 Strip, 0.2 mL			1 Low Veltage Double Insulated Down	
Sprout [®] Plus Lid O-Ring, Pack of T	10		 1 Low Voltage, Double Insulated Power 240 VAC, 50/60 Hz) with 4 Interchange 	
Sprout [®] Plus Power Adapter				
Tube Adapters, 0.2 mL, 0.5 mL, P	ack of 12 (6 ea.)			

Sprout[®] Plus Power Adapter Tube Adapters, 0.2 mL, 0.5 mL, Pack of 12 (6 ea.)

Specially Designed Inserts seat the centrifuge, all its accessories, and includes a tube rack



Portable Centrifuge Kit

Spin samples on-the-go

The Portable Centrifuge Kit offers a convenient way to spin down and prepare samples outside of a conventional laboratory environment

Immediate Field Processing

- Process samples in remote areas before deterioration can take place
- Allows for immediate on-site testing for faster results
- Made of durable ABS plastic for easy clean in field environments

Easy to Use

- Centrifuge is equipped with tool-less SnapSpin rotors for easy change out
- Simple operation with lock/unlock activation

Item No.

120568	Portable Centrifugation Kit, 100-240 VAC, 50/60Hz Universal Plug + 12V Car Adapter
120569	Portable Centrifuge Transport Case with Inserts - Only
120570	12V Sprout [®] Plus Centrifuge Car Adapter*
	*Adapter can be used with Sprout® Plus Mini Centrifuge

Portable Convenience

- All items in kit are housed in a handled DuraPorter[®] Transport Box for portability
- The handled box has an air-tight seal for optimal sample protection



12-Volt Car Adapter

for on-the-go testing

Package Includes

Sprout [®] Plus Mini-Centrifuge
6-Place 1.5/2.0 mL Tube Rotor
16-Place 0.2 mL Tube Rotor (singles, or 2 x 8 strip)
6 Tube Adapters 0.2 mL
6 Tube Adapters 0.5 mL
4 Lid O-Rings
Storage Case for Rotors and Adapters
Tube Rack
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs
DuraPorter® Transport Box
Tube Rack Insert
(12-Place 1.5/2.0 mL, 14-Place 0.2 mL)
Transport Box Dividers

Sprout[®] Plus Cradle

Portable 12V Sprout® Plus Centrifuge Car Adapter

Magnetic Induction Stirrer

Strong, motor-less and maintenance free

Superior Functional Design

- Advanced magnetic coil technology drives the spinning of the stir bar maintaining continuous, constant speed for reproducible results
- Magnetic coupling is optimized by using alternating current (AC) as its driving force strong enough to stir viscous solutions
- Motor-less operation generates minimal heat ideal for stirring temperature sensitive samples and environments for long periods
- With no moving parts, produces less vibration and quieter operation over traditional magnetic stirrers

Intuitive Operation

- Multiple stirring modes are easy to program
- Bright digital display precisely controls time and speed
- Quickly meets stirring speeds with rapid acceleration and braking (less than 10 seconds)



Item No.	Color	L x W x H in	L x W x H cm
120584	Grey/Purple	9.7 x 7.7 x 3 24.5 x 19.5 x 7.5	04 E x 10 E x 7 E
120585	Grey/Blue		24.3 X 19.3 X 7.3

Replacement	parts
120738	MI Stirrer Silicone Mat
120739	MI Stirrer Universal Power Adapter



Protective silicone anti-skid plate cover for improved safety





Durable ABS housing is corrosion resistant and easy to clean

Specifications Speed Range 50–2,000 rpm variable, 1 rpm increments Run Time Range 30 seconds to 60 minutes, variable, 1 second increments or continuous Stirrer Modes Single clockwise Single counter clockwise Single counter clockwise

	-			
	Auto reverse rotation	1 to 99 cycles 50% duty		
Capacity	3 L plus			
Mixing Surface	165 mm diameter, silio	165 mm diameter, silicone mat 160 mm diameter		



Package Includes	
Rare Earth, NdFeB Magnetic Induction Stirrer	
2 Stir Bars 35 mm x 12 mm Cylindrical Shape	9
1 Silicone Anti-Skid Plate Cover	
1 Low Voltage, Double Insulated Power Adap Interchangeable Plugs	ter (100-240 VAC, 50/60 Hz) with 4

9



Mini Magnetic Stirrer

Affordable solution, ultimate function

The Mini Magnetic Stirrer is an affordable, dependable solution ideal for general research, academic, life sciences, industrial and clinical lab settings

Durable Design

- Constructed of chemical resistant ABS and PC blend
- Compact footprint saves valuable space on benchtop
- Lightweight for easy transport

Ultimate Function

- Powerful magnetic coupling ensures accurate alignment, reducing spin-outs during acceleration or deceleration
- Performs routine stirring applications with low vibration and quiet operation



Specifications

U.S.and E.U. Utility Patent Pending

Speed Range 200–2,200 rpm (Infinitely adjustable)

Capacity Up to 3 L

Package Includes Mini Magnetic Stirrer



			U.S. Design Patent No. D814,045
Item No.	Color	L x W x H in	L x W x H cm
120155	Grey/Purple 🔳 🔳	7.2 x 5.2 x 1.85	18.3 x 13.2 x 4.7
120595	Grey/Blue	7.2 X 3.2 X 1.03	10.3 X 13.2 X 4.7


Replacement Pa	rt
120815	Mini Stirrer Universal Power Adapter



2 Stir Bars 35 mm x 9 mm, Octagonal Shape 1 Low Voltage, Double Insulated Power Adapter

Stir Bars & Stir Bar Retriever

- Magnetic bar is made from ALNICO, one of the strongest permanent magnets available
- The sintered manufacturing process reduces the chance of decoupling (spin outs), encapsulation in PTFE protects both magnet and what it is stirring
- Chemically inert and easy to clean
- Can withstand temperatures up to 250°C, making them suitable for use with both unheated and heated stirrers
- Low friction octagon shape with a center pivot ring helps create a strong mixing vortex, even with more viscous solutions
- The retrieval bar can recover stir bars even in hot solutions when stirring and mixing is no longer required

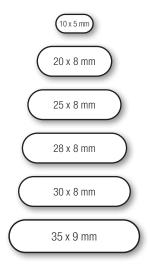

Stir Bars		Pack of 10
ltem No.	L x W (in)	L x W (<i>mm</i>)
120547	0.4 x 0.2	10 x 5
120548	0.8 x 0.3	20 x 8
120549	1 x 0.3	25 x 8
120550	1.1 x 0.3	28 x 8
120551	1.2 x 0.3	30 x 8
120552	1.4 x 0.4	35 x 9

Stir Bar Retri	Pack of 5	
Item No.	L (in)	L (cm)
120553	11.8	30







Vortexer

Multi-purpose, maximum convenience

Innovative Design

- One mixing head does it all with this sleek, modern and compact designed vortex mixer
- Hold slot configurations for multiple sizes of tubes, vessels and assay plates
- Center pad for quick hand-held mixing

Easy to Use

- Control switch with three operation selections (intermittent, off, continuous)
- Variable Speed Rotary Dial provides pre-set speed control selection during shaking mode

Durable Construction

- ABS material resist most commonly used chemicals and is easy to clean
- Stable elastomer feet dampen vibration and _ prevent movement during operation




hand held mixing



5 mL tubes x 2 (horizontally mounted in the same retention space as 15 mL)



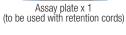



0.5 mL microfuge tubes x 6 (vertically mounted)



25 mL, 50 mL conical tubes x 2 (horizontally mounted)

(vertically mounted)



(horizontally mounted)



0.2 mL microfuge tubes (strip tubes) x 24 (vertically mounted)

US Design Patent No. D787,085 US Utility Patent No. 9895670

Item No.	Speed Control	Power Input		L x W x H in	L x W x H cm
120209	Infinity Dial	110-120 VAC, 50/60 Hz, 55 W, US/NA	1,1		
120212	Infinity Dial	220-240 VAC, 50 Hz, 44 W, CA	$\langle \cdot \cdot \rangle$	6.7 x 7.5 x 7.8	17.1 x 19 x 19.7
120214	Infinity Dial	220-240 VAC, 50 Hz, 44 W, UK	<u>.</u>	(with mixing head	(with mixing head
120210	Preset Dial	220-240 VAC, 50 Hz, 44 W, AUS	\odot	in place)	in place)
120318	Preset Dial	100-110 VAC, 50/60 Hz, 66 W(PSE), Japan		-	

Replacem	ient parts	Replace	Replacement parts		
120327	Vortexer Multi-Head	120333	220-240 VAC AUS Cord	\odot	
120328	Vortexer Central Pad, 3 pk	120332	220-240 VAC CE Cord	$\langle \cdot \cdot \rangle$	
120422	Vortexer Retention Cords, 3 pk	120331	220-240 VAC UK Cord	<u> </u>	
120329	100-120 VAC US Cord	120330	100-110 VAC Japan Cord		

Specifications:			
Speed Range Approx. 1,000; 2,000; 3,000 rpm			
Orbit	3.7 mm		
Control	3 Position Control — (top) intermittent, (middle) off, (bottom) continuous		



Models 120209, 120212, 120214



Mini Vortexer Small footprint with maximum power

Simple, touch start operation offers exceptional mixing performance with one hand

Robust Design

- Low profile design and small footprint saves valuable space
- Made of ABS and silicone materials that resist most commonly used chemicals
- Stable platform and elastomer feet ensure unit stays stable during operation

Powerful Function

Replacement parts

120617

120814

- Rapid vortex action with low activation force
- Quiet operation and low vibration





Item No. Color LxWxH LxWxH 120567 120598 120827

Speed	4,500 rpm
Orbit	6 mm
Capacity	Up to 50 mL tube

		<u>III</u>	6111
7	Grey/Purple		
3	Grey/Blue	5.1 x 4.3 x 3	13 x 11 x 7.5
7	Black/Black		

	• · · · · • • • •	
Black/Black		
		 Package Includes
ent parts		Mini Vortexer
Mini Vortexer Silicone	Cover, set of 2	2 Silicone Caps
Mini Vortexer Universa	I Power Adapter	

1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs

$(\bullet \bullet)$	 (\bigcirc



Digital Orbital Shaker

Shaker for a vast range of protocols

Versatile Design

- Create your own matrix for vessels of all shapes and sizes with Hi-Lo Platform and Bungee Cords
- Instantly adjusts to weight load changes and unbalanced loads with nostop electronic motor feedback allowing for truly reproducible results
- Durable ABS construction suitable for use in incubators (excluding C02) and cold rooms within the stated temperature range

High-Featured Operation

- Suitable for wide range of protocols including very low start speed for sensitive sample agitation to high speeds for vigorous applications
- Soft start and shut down feature reduces spilling from full vessels
- Designed to handle heavy loads up to 10 lbs (4.5 kg)
- Overload protection shut-off prevents overloading

Item No.	L x W x H in	L x W x H cm	Platform Int	terior W x D cm
120460	12.4 x 7.8 x 12.2	31.4 x 19.7 x 31	11 x 11	27.9 x 27.9

Replacem	ient parts
120521	Digital Orbital Shaker Silicone Mat & Bumpers
120524	Digital Orbital Shaker Bungie Cords, Pack of 8
120491	Digital Orbital Shaker Universal Power Adapter



Specifications	
Speed Range	20–300 rpm, 1 rpm increments
Orbit Diameter	19 mm (0.8")
Run Time Range	0 to 74 hours and 59 minutes; 1 second increments or continuous
Maximum Load	4.5 kg (10 lbs)



Package Includes
Digital Orbital Shaker
Skid-Free Rubber Platform Mat
8 Bungee Cord Set
Universal Hi-Lo Platform with Soft Rubber Edge Bumpers
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs Interchangeable Plugs



Easy Identification

5 different colored nose cones help differentiate department, task or user; as well as a labeling area on the side of the unit



Magnetic & Screw Wall-Mount



3 Modes Bright LEDs indicate 3 operation modes: high, low, & gravity



Specifications	
Power	Micro-USB connector Rechargeable li-Ion battery (replaceable)
Filter	PTFE/25mm dia/0.45 µm filter
Pipettes	1–100 mL, glass and plastic
Pipette Nozzle	Polypropylene



RF3000 [®] Pipette Controller		
Charging Stand and Wall-Mountable/Magnetic Holder		
Nozzle Adapter for use with Small Pipettes		
Li-Ion Rechargeable Battery		
Five Nose Cones: Black, Yellow, Red, Grey and Blue		
5VDC Micro USB Plug		
6 ft (183 cm) Cable		
1 Low Voltage, Double Insulated Power Adapter (100-240 VAC, 50/60 Hz) with 4 Interchangeable Plugs II ••• I		

RF3000®

Packed with features for comfortable and accurate pipetting

The powerful pump offers rapid, yet sensitive aspiration and dispensing, has adjustable valves and is microprocessor controlled for constant reliability and accuracy

- Long-life li-ion battery, with quick charge-up to 20 hours use
- Toggle operation by pressing the large mode button
- An adapter can be fitted to the nose cone to stabilize small volume and thinner pipettes



the unit

Nose cone assembly only

- Low-battery indicator

- New USB connecter ensures the

- Can be recharged when in use

appropriate charger is used with



doubles as a charging station

Item No.	Volume Range mL	L x W x H in	L x W x H cm	
HS3000	1–100	9 x 8.5 x 5	22.9 x 21.6 x 12.7	
Replacement p	arts			
HS5000RFFIL	RF3000 [®] Replacement Filters, <i>Pack of 5</i>			
120669	RF3000 [®] Universal Power Adaptor 5VDC Micro USB Plug			
HS300502	RF3000® Li-Ion Rechargeable Battery, One year warranty on battery included			
120487	RF3000® Nose Cones 🛛 Assorted* 🔳 🔳 🔳 📕	Set of 5		
HS5000RFRFA	RF3000 [®] Replacement Filter Adapter			
HS5000RFRP	RF3000 [®] Replacement Pipette Adapter			
		*Assorted colors (Black, Blue	, Grey, Red and Yellow)	

Pipette Controller

Precise fluid control at your finger tips

Comfortable to hold and easy to use, the manual inline pipette controller was designed with an ergonomic grip and bulb for simple and precise operation

A T-bar thumb lever controls the aspirating and dispensing speeds from serological pipettes from 1 to 100 mL in size

- Conveniently accommodates either glass or plastic serological or volumetric pipettes ranging from 1–100 mL
- Includes a replaceable .45 µm hydrophobic filter to protect your samples and internal components



Nose cone assembly only

Item No.	Color	Volume Range mL	L x W x H in	L x W x H cm	
120479	Blue 🗖	1–100	9 x 2.1 x 2.1	23 x 5.4 x 5.4	

Operation Is Simple:

Squeeze the bulb to create a vacuum to draw up fluid

Push the control lever forward (down) to dispense the liquid

Push the control lever backward (up) to aspirate liquid through a pipette

Press the central bulb on the unit to blow out any remaining fluid if necessary



Stay Organized with two additional nose cones (green, white) for color identification for both the Pipette Controller & RF1000



RF1000

Lightweight and easy to use

The RF1000 is a manual pipette controller that is uniquely designed for one-handed operation

Simply squeeze the main bulb and push the toggle forward to aspirate or back to dispense. Push the blowout bulb to expel any liquid remaining in the pipette

- Tapered silicone rubber pipette adapter allows a large range of pipette sizes ranging 1-100 mL, replaceable when worn
- Includes a replaceable .45 µm hydrophobic filter to protect your samples and internal components



Nose cone assembly only


Item No.	Color	Volume Range mL	L x W x H in	L x W x H cm
HS1000RFU	Blue 🗖	1–100	7.1 x 1.8 x 4.7	18 x 4.5 x 12

Replacement Parts

Item No. Description		Color
HS5000RFFIL	Replacement Filters, Pack of 5	
120604	Nose Cone with Adapter	Blue 🗖

Squeeze to create a vacuum TITITITITITITITITITITITI



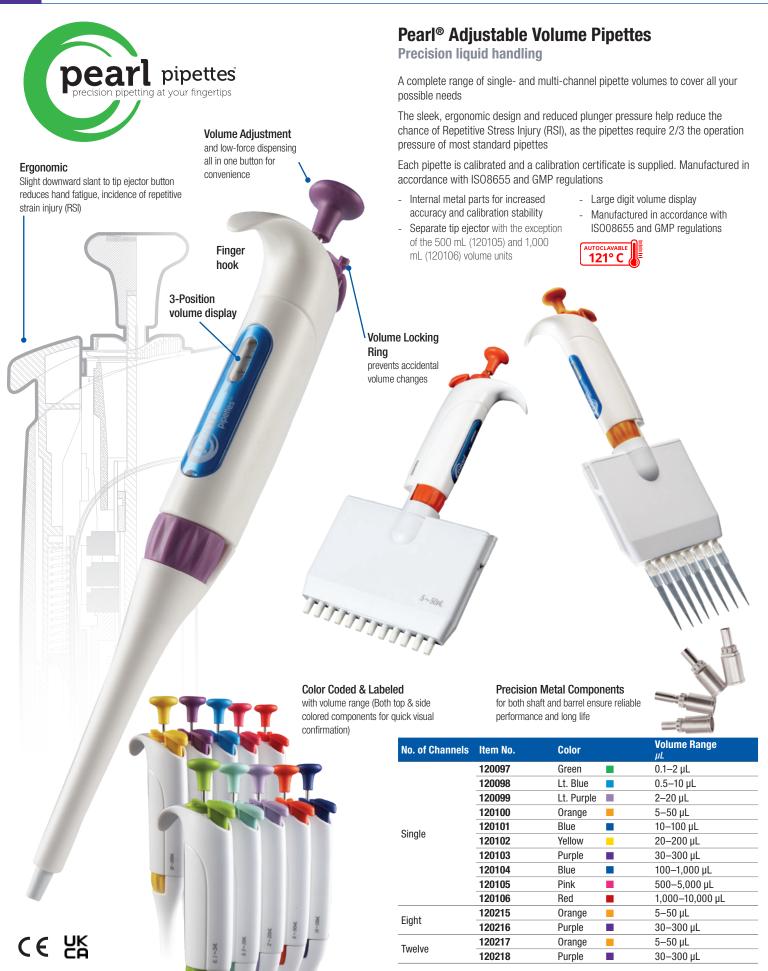

Pipette Starter Kit

Universal Carousel Pipette Stand

Hold and store pipettes in style with the Universal Carousel Pipette Stand by Heathrow Scientific that holds up to 8 single-channel, 4 multi- channel, or a combination of both pipettes

A sleek and modern design enables this stand to securely hold virtually every known brand of pipette in the marketplace making it truly 'universal' while adding a pop of color to a bench top

Reagent Reservoirs

Multi-channel reservoirs by Heathrow Scientific are superior in every way with thicker walls, crush-resistant "turrets" at the corners, easy pour-spouts, 25 mL graduation marks, and tabs to prevent stacked reservoirs from sticking together

Pearl® Pipettes

Pearl Pipettes by Heathrow Scientific were specifically engineered for ease of use and comfort

The sleek, ergonomic design and reduced plunger pressure help reduce the chance of Repetitive Stress Injury (RSI), as the pipettes require 2/3 the operation pressure of most standard pipettes

Each pipette calibrated and a calibration certificate is supplied

Manufactured in accordance with ISO8655 and GMP regulations

Pipette Starter Kit Item No. 120554

Universal Carousel Pipette Stand

Upper Holder Color	Lower Holder Color	Base Width in	Top W x H in	Base Width cm	Top W x H cm
Assorted/Clear	Clear 🗌	6.2	9 x 16.1	15.8	22.9 x 40.8

*Assorted colors (Purple, Blue, Green, Grey, Red Orange, Yellow and Clear)

Pearl® Pipettes - Single Channel

Color			Volume Range µL				
Lt. Purple			2–20 µL				
Yellow			20–200 µL				
Blue			100–1,000 µL				
Reagent Re	eservoirs						
Color		Volume mL	Material	Sterility	L x W x H in	L x W x H cm	
White		50	Polystyrene	Sterile	4 x 1.1 x 1.2	10.2 x 2.8 x 3	Pack of 5

*All parts available separately





Rapid, accurate colony counter

Combines a multi-function electronic counter with a Sharpie[®] Pen for marking to prevent missing or double-counting colonies

- Light, ergonomic body is balanced for easy manipulation
- Select between an audible beep or an LED flash to acknowledge count
- Store up to 32 separate counts and display the total of all counts

Package Includes:

- eCount Colony Counter
- Black Sharpie® Pen
- Long-life button battery included

- Count up/down verification
- Quick and easy to change
 - between different colored markers to check different colony counts in the same view

- Hand tally AdapterWeighted Stand

US Patent No. D564,385 EP Registered Design No. 705751-0001

	-	
Item No.	L x W	L x W
	in	ст
120000	7.4 x 1.2*	17 x 3*
120228	CR1225 Battery	

*Measurements shown are with Sharpie® Pen inserted






Weighted Stand for bench top use







Lock Counter Off to allow use of the marker without registering more counts

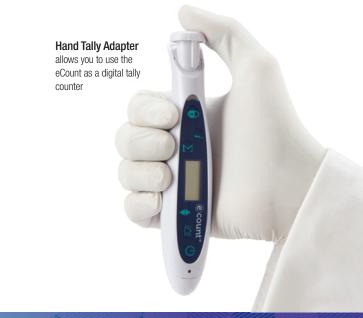



Mini-Light Box and Mini-Magnifier

Perfect accessories for the eCount or use as stand-alone items

Mini-Light Box

- LED Light source. Color temperature 5,300K +/- 270K
- Viewing area (125 x100 mm) is perfect for quick screening of 96-well
- plates, petri dishes, electrophoresis gels, slides, and transparencies
- Powered by six 1.5V AAA batteries (included), not included in warranty

Mini-Magnifier

- Swiveling molded lens 1.75x magnification
- Enlarges objects for easier viewing

Item No.	Description	L x W x H in	L x W x H cm
120661	Mini-Light Box*	6.2 x 5.9 x 0.6	15.8 x 15 x 1.4
HS15092D	Mini-Magnifier	6.3 x 5.3 x 8.1	16.1 x 13.5 x 20.7
			*Batteries included



Hand Tally Counter

- Features a quick-reset knob

Use for counting blood cells, bacterial colonies, drops of liquid or any repetitive event

- Tallies up to 9999
- Lens for easy reading, and finger ring for a secure grip
- Item No.
 Dia. x W in
 Dia. x W cm

 HS6594
 1.8 x 1.6
 4.6 x 4.1



Benchtop Counters

Mechanical, benchtop counters are perfect for blood cells (label included), bacterial colonies or repetitive counting

- Each key records up to 999
- Blood cell key labels
- Total window keeps track of total number of strokes made on all other keys; bell rings when total reaches 100

Item No.	No. of Keys	L x W x H in	L x W x H cm
HS5360A	6-Key	9.8 x 3 x 2.2	25 x 7.6 x 5.6
HS5360B	9-Key	12.5 x 3 x 2.2	32 x 7.6 x 5.6

